

LARYNGOSCOPES

HANDLES, CONVENTIONAL CHROME-PLATED BRASS

A complete range of conventional (bulb in blade) laryngoscope handles expertly fabricated from brass, polished, and durably plated with chrome. Handles are available in a variety of sizes and configurations to assist in both routine and extremely difficult intubations.

STANDARD (29mm x 140mm)
Lightweight standard. Requires two "C" size batteries (#P-154-C). **#00-350**
Standard size in sealed/waterproof configuration. **#00-350-WP**

PENLITE (19mm x 140mm)
Slender profile. Requires two "AA" size batteries (#P-154-AA). **#00-354**
Penlite size in sealed/waterproof configuration. **#00-354-WP**

PEDIATRIC (16mm x 100mm) **#00-336**
Tiny and short. Requires two "N" size batteries (#P-154-N).

LARGE (38mm x 170mm) **#00-341**
Large and stout. Requires two "D" size batteries (#P-154-D).

BANTAM (32mm x 100mm) **#00-342**
Shorter and lighter than the Standard. The Bantam handle is designed to accommodate the difficult intubation on a short neck, barrel-chested or large bodied individual where the distance between the mouth and the chest prohibits the use of a standard length handle. The size is such that it accommodates the normal sized hand and can be used on a routine basis by most people. Requires two "AA" size batteries (#P-154-AA).

ADJUSTABLE (32mm x 130mm) (Patil-Syracuse) **#00-340**
Similar in size and weight to the Standard, yet readily allows the blade to be firmly locked at four different angles for a variety of special uses:
45° - Howland Lock position, natural lift for bull-neck, anterior larynx
90° - Standard position, identical to a Standard handle
135° - Polio blade position, clears obstructions, extra room
180° - Parallel to handle, for obese, halo-traction, or cervical fractures
Requires two "AA" size batteries (#P-154-AA).

ACCESSORIES, CONVENTIONAL

HOWLAND LOCK **#00-343**
Natural lifting action simplifies laryngeal exposure, making intubation possible even in the most difficult cases (receding chin, anterior larynx, protruding teeth, "bull neck", facial constrictures, decreased jaw mobility, etc.). Built-in leverage prevents prying and reduces possibility of broken teeth. Adapts to all types of conventional laryngoscope handles and blades. Can often make an otherwise certain tracheotomy unnecessary.

REPLACEMENT LAMPS


Standard (frosted) pkg of 12	Large	#00-347
	Small	#00-348
Pre-focused (clear) pkg of 5	Large	#00-355
	Small	#00-356

Extra Bright (clear) pkg of 5	Large	#00-355-XB
Superior gas charge in bulb, parabolic reflector on large size, 400% more light output relative to standard	Small	#00-356-XB


REPAIR PART / KITS

Blade rewire kit w/bulb, present AincA and similar	Large	#00-368
	Small	#00-369
Blade rewire kit w/bulb, Penlon style	Large	#00-368
	Small	#00-369
Blade rewire kit w/bulb, old style AincA	Large	#00-338
	Small	#00-339
Handle contact assembly, 3/4" length		#00-371
Handle pivot pin, screw type		#00-372

Note: Other types of Specialty Blades may be available on special order.

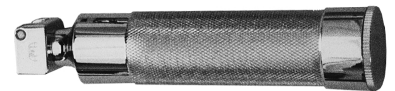


#00-350 Standard

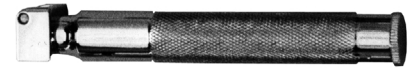


#00-354 Penlite



#00-336 Pediatric

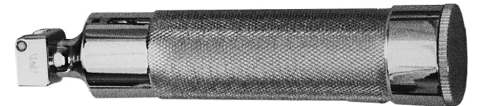


#00-341 Large

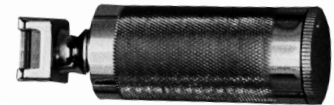



#00-342 Bantam



#00-340 Adjustable

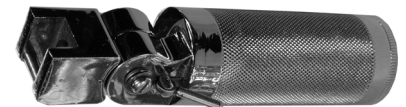

#00-343 Howland Lock

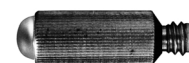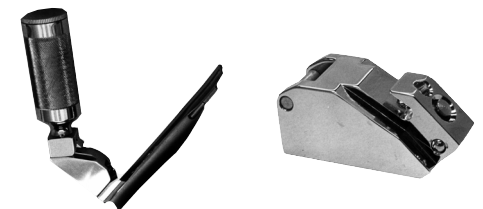

Large Lamp



Small Lamp