

# PRESSURE MANOMETERS AND ACCESSORIES

**Manometer #00-199**  
 Manometer has range of -40 to +80 cm H<sub>2</sub>O in a 2-1/2" diameter chrome-plated case with 1/4 MPT bottom connection. Has front zero correcting screw.

**Manometer, Max Hold Pointer #00-199-MAX**  
 #00-199 gauge with resetable low friction pointer to capture maximum pressure or vacuum reading.

**Manometer, MRI Compatible #00-199-MRI**  
 Similar to #00-199, but constructed of brass and synthetics to allow use in powerful MRI environment. No ferrous components.

**Gauge Guard (only) #00-150**  
 To fit #00-199 style Manometers. Constructed of PVC elastomer with rear mounted attachment strips.

**Dovetail Manometer Mounting #00-195**  
 For use with #00-199 or #00-274 Manometer. Male portion of dove-tail mounts to any flat surface by means of screws provided. Female portion attaches to mounting screws on back of gauge.

**Swivel Manometer Mounting #00-196**  
 For use with #00-199 or #00-274 Manometer. Mounts to any flat surface and to the base of the gauge. Provides 4 1/2" swivel arms w / 180° dual pivot. Stainless and nylon.

**Pole Clamp Manometer Mounting #00-200**  
 For use with #00-199 or #00-274 Manometer. Accommodates up to 1" pole, square or round. Attaches to mounting screws on back of gauge.

**Pressure Relief Valve (APL) #00-273**  
 True APL adjustable valve with settings of 0-55 cm H<sub>2</sub>O at 9 lpm. Mounting thread is 1/2-20.

**Manometer w / APL Valve #00-274**  
 #00-199 Manometer above with P-159 "Tee" (1/4" tube) and #00-273 APL valve adjustable from 0-55 cm H<sub>2</sub>O at 9 lpm.

**Manometer, Max Hold w / APL Valve #00-274-MAX**  
 #00-274 with resetable pointer to capture maximum pressure or vacuum reading. (Not shown)

**Manometer w / Harness and Adapters #00-275**  
 #00-199 Manometer above with 15mm female adapter #00-275-1, neoprene elbow #00-452, and harness #00-276. Used post-operatively for lung function evaluation.


**Manometer, Max Hold w/ Harness #00-275-MAX**  
 #00-275 with resetable pointer to capture maximum pressure or vacuum reading. (Not shown)

**Adapter, Tube #00-489**  
 Threaded tubing adapter for use with #00-199 Manometer. 1/4 F.P.T. by 1/4 O.D. tubing attachment.

**Adapter, Tee #P-159**  
 "Tee" 1/4 F.P.T. x 1/2-20 with 1/4 O.D. tubing attachment.

**Replacement Crystal, Standard #P-158**  
 For #00-199, #00-274, and #00-275 gauges.

**Replacement Crystal, Max Hold Pointer #P-158-MAX**  
 For #00-199-MAX, #00-274-MAX, and #00-275-MAX gauges.

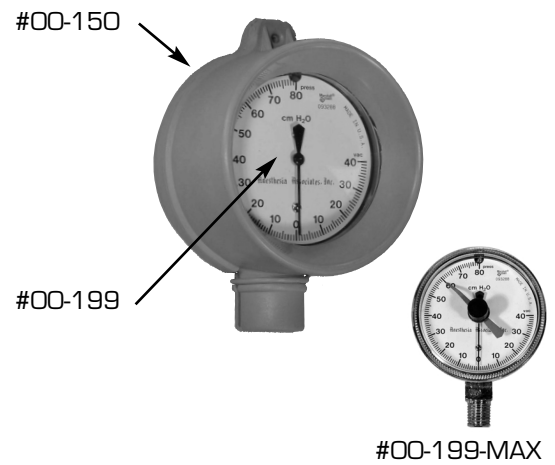

#00-199-MAX



#00-195

#00-196

#00-200



#00-274



#P-159




#00-273



#00-489

#00-275

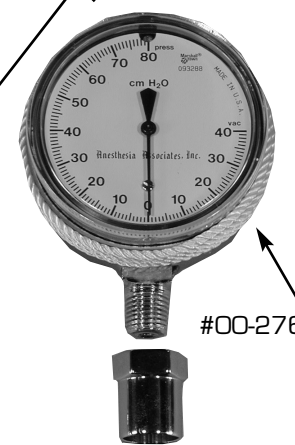

#00-276



#00-452

#00-275-1